




# WHITE PAPERS & APPLICATION NOTES

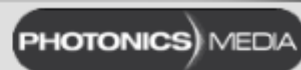

**DOWNLOAD FREE WHITE PAPERS & APPLICATION NOTES**

## Working in the Basement: Measuring Signals Below the Noise Floor with a Lock-In Amplifier

Measuring optical signals in the femtowatt (10<sup>-15</sup>) to nanowatt (10<sup>-9</sup>) range can be a daunting task. Signal levels this low are lost in typical detector noise levels and swamped by background light. The noise floor for photodiode detectors operated with a small bandwidth (~10 Hz) is on the order of 1 picowatt (10<sup>-12</sup>). In order to achieve significant improvements in noise rejection we need to turn to a lock-in amplifier. Lock-in amplifiers can improve noise rejection by 3 orders of magnitude or more.

[DOWNLOAD NOW](#)

**Working in the Basement: Measuring Signals Below the Noise Floor with a Lock-In Amplifier**  
By Shimon Elstein, Senior Physicist, Ophir

Measuring optical signals in the femtowatt (10<sup>-15</sup>) to nanowatt (10<sup>-9</sup>) range can be a daunting task. Signal levels this low are lost in typical detector noise levels and swamped by background light. The noise floor for photodiode detectors operated with a small bandwidth (~10 Hz) is on the order of 1 picowatt (10<sup>-12</sup>). Further narrowing of the bandwidth by filtering or averaging will only provide a small additional reduction in the noise level.

In order to achieve significant improvements in noise rejection we need to turn to a lock-in amplifier. Lock-in amplifiers can improve noise rejection by 3 orders of magnitude or more. Furthermore, they can provide background signal rejection that is several orders of magnitude higher than the noise rejection.

Lock-in amplifiers employ a technique called homodyne detection to obtain their outstanding performance. Homodyne detection has two requirements: a) the signal to be detected needs to be modulated and b) a clean reference signal with the same frequency needs to be provided. In the lock-in amplifier, the signal to be measured is multiplied by the reference signal and then integrated over time. This results in an extraordinarily narrow effective bandwidth. Signals at all frequencies that differ by even a small amount from the reference frequency will result in a net integrated value of zero. Since detector noise is "white", its power is spread across a broad spectrum and the noise amplitude component at the measurement frequency is very low. By confining the measurement to a single frequency, detector noise is reduced drastically. In a similar fashion, background optical signals (primarily DC or line frequency) are rejected by the lock-in amplifier.

The key to high performance with a lock-in amplifier is maintaining a precise match between the modulation frequency of the signal to be measured and the frequency of the reference signal. In optical applications this precise match is readily built in to the architecture. Many low level optical signals that need the benefits of a lock-in amplifier are DC or very low frequency. In these applications, an optical chopper is used to modulate the signal.

Figure 1. Typical Optical Chopper Disk

Sponsored by



### More White Papers from this Sponsor

- Reducing Production Waste with Laser Profiling and Characterization
- Laser Measurement Best Practices: How and When to Apply Measurements
- Measuring Non-Light Sources with Commercial Laser Beam Profiling Systems

## PHOTONICS MEDIA

Visit Photonics Media to download other white papers and learn more about the latest developments in lasers, imaging, optics, biophotonics, machine vision, spectroscopy, microscopy, photovoltaics and more.

[www.photonics.com/WhitePapers.aspx](http://www.photonics.com/WhitePapers.aspx)

Questions: [info@photonics.com](mailto:info@photonics.com)

[Unsubscribe](#) | [Subscribe](#) | [Subscriptions](#) | [Privacy Policy](#) | [Terms and Conditions of Use](#)

Photonics Media, 100 West St., PO Box 4949, Pittsfield, MA 01202-4949  
© 1996 - 2017 Laurin Publishing. All rights reserved. Photonics.com is Registered with the U.S. Patent & Trademark Office.  
Reproduction in whole or in part without permission is prohibited.